TIANYI-BIOM
MICRO-BIOENERGY VITALIZER

User Manual

Product Description

Building on the foundation of Pulsed Electro Magnetic Field (PEMF), the Micro-Bioenergy Vitalizer is a successful breakthrough as a result of years of research and a combine effort from professionals in different fields, ranging from elite health, biomedical science to electronics engineering and information technology.

As it advances to the third generation, it performs many functions. Not only does it stimulate cellular and bioenergy regeneration, do pulse-check, it also performs as massager. By sending signals to cell receptors to regenerate, combine with the vibration on the meridian points, it combines the traditional TCM and technology to improve blood circulation and boost organ functions, resulting in significant health benefits.

With this non-invasive, easy-to-use technology, and constant effort from the user-end, improving health is just a few massages away! The ultimate aim is to help people improve their blood circulation, bioenergy, and feel revitalized anywhere not limited by time and space so it can be accessible anytime, anywhere.

1. Safety Instructions

1. This product has passed EMC and electrical safety test for CE market.
2. Power: built-in 3.7V lithium battery power supply, external charger;
3. Input voltage: input DC5V $\overline{\text{---}}$ 1A;
4. Battery capacity: 1500mAh;
5. Overall size specification: elliptical column: 150 × 39 × 29mm, \pm 2mm;
6. Pulse rate display range: 30--250bpm, \pm 2bpm;

2. Functions

The Micro-Bioenergy Vitalizer as two main functions, pulse checking and meridian points treatment.


1. Pulse checking: place left middle finger on the sensor behind, results shown in the display represents pulse per minute.
2. Massage treatment: place device on positions as directed on the display.

3. Expected Outcome

- Blood-circulation, and improved organ functions.
- Re-energize, revitalize cells;
- Blood cleanse, improved wellness.

4. Pictorial Description

As per below, main body consists of display window; power button; start/stop/pause button; signal point; and sensor.

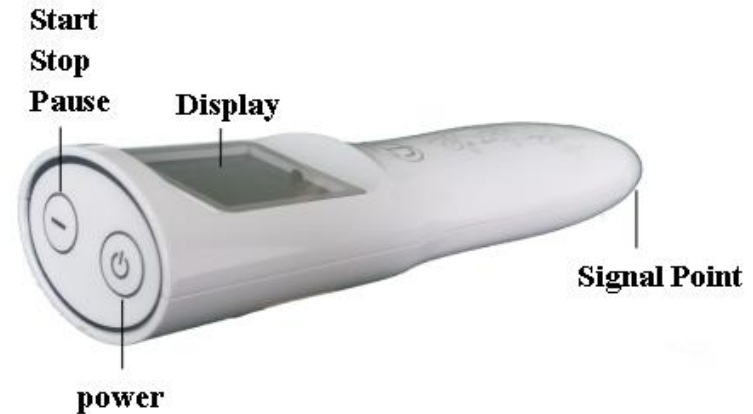power

Pulse Sensor

5. Meridian Positions

Place device as per pictures shown below and maintain the positions until timer stops.

Mid-palm 5 min; 3-finger from wrist 3 min.

6. Step-By-Step Guide

2 modes:

1. Pulse checking;
2. Direct treatment mode.

Details:

1. Pulse Checking

1. Place left middle finger on the sensor located at the back, press the power button for 3 sec. Screen will turn on as per picture shown below.

2. Display will show time and testing page.

3. Pulse rate will be obtained within 30 seconds and display as per picture below. The pulse shown is 80bpm.

4. In treatment mode, align the signal point to the middle of palm, press the start button. Hold position for 5 minutes.

5. Change position upon seeing the following screen.

6. Place the device 3 fingers spacing down from wrist and press start button. You should see the following:

7. Upon timer reached, the treatment is complete. The instrument will auto-shutdown in 3 seconds thereafter.

2. Direct Treatment:

1. Place device in the mid-palm position and press power button, it will show start up screen as shown below. Press start button to skip pulse checking.

2. Upon pressing start button, the screen will show as per picture below.

3. Change position upon seeing the following screen.

4. Place the device 3 fingers spacing down from wrist and press start button. You should see the following:

5. Upon timer reached, the treatment is complete. The instrument will auto-shutdown in 3 seconds thereafter.

7. Important

1. Do not use the device on the same spot within 2 hours after use;
2. This product is a massager for meridians and collaterals, and does not replace any medical treatment. Please do not use it if a pacemaker has been installed;
3. The 5th use of the day should be before 1130pm;
4. Do not place device near flammables;
5. Never use the device under extreme hot environment;
6. Do not turn the device on/off repeatedly, keep it away from children;
7. Do not use device under drunk or extreme fatigue conditions, use after rest;
8. If the device dropped, please bring to service center to check its functionality;
9. Never try to disassemble the device;
10. Do not use when charging;
11. Do not share personal-configured device with others;
12. Any queries, please contact customer service;
13. Pulse readings do not interfere the treatment;
14. **If needed, press the pause and resume the treatment by pressing it again.**
15. To shut down, press the power button for 3 seconds.

16. During treatment, the following screens should display.

8. Charging

Connect the charger to 220V power supply, look out for hazard, then connect the USB cable to device. Blue LED indicates charging. Per charge is 2-3 hours, do not prolong charging time.

9. Other Information

1、Display details

2、Low battery

(The following displays when battery falls below 20%)

3、Pulse Reading Error

(If pulse reading is unable to obtain within 30 sec, as shown below, please restart the device and try again)

Heart rate can not de detected, try again or press the start button directly.

Heart rate can not de detected, try again or press the start button directly.

4、System Error

(If encounter system error, as shown below, please restart the device and try again. If the message persists, contact customer service for assistance)

5、Overuse

(Please wait 120 minutes before you can use it again)

10. Support

Warranty

Dear user:

Thank you for your support, the following services are available to you upon purchase:

1. Under normal operation, one year hardware coverage and life-time technical support.
2. Any changes or modification to the product would result in warranty void and repair charges apply.
3. Please bring warranty card and receipt for repair.

Warranty Coverage

1. One-year product hardware coverage,
2. Courier charges shall be paid by company for repair;
3. Charges apply for the following:
 - a. Improper use of the device;
 - b. Device disassembly, out of warranty;
 - c. Other damage caused by environment or users;
 - d. Mishandling of equipment by user or modification;
4. Present warranty card for repair.

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Copyright
Tianyi Biomedical Tech Pte. Ltd. Singapore